



ENMED® ENM - 224 STERILE with Towel & ENM - 214 NON-STERILE (ENM - 2x4) with Towel

PRODUCT DESCRIPTION: The ENMED™ reinforced surgical gown (ENM-224 sterile and ENM-214 non-sterile) is manufactured from a breathable, nonwoven, polypropylene fabric, with SMS material at the back. ENM-224 is produced with reinforced SSMS fabric on the chest and arms, which are critical areas of the gown. The film layer is a multilayer polypropylene olefin, breathable film (SSMS). The ENMED™ reinforced surgical gown variants have been tested according to AAMI PB70:2012 and meet LEVEL 2/3/4 requirements of the AAMI liquid barrier classifications for a surgical gown.

The ENMED™ reinforced surgical gown (ENM-224 sterile and ENM-214 non-sterile) is closed-back and made out of laminate with ultrasonic seams. It has a velcro closure at the back of the neck and a waist tie feature to secure the gown to the body of the user.

There are four belts (two inside two outside the gown) made of the same material as the gown and ultrasonically welded to the gown. The belts are used to secure the closure of the back.

A ribbon tie from each side is provided as extra support to close the neck.

The sleeves of the gown have knitted cuffs (made out of polyester) sewn onto the end of the sleeves at the user's wrists to keep the sleeves in place on the wearer.

The entire gown, except for the ribbon tie and knitted cuffs, are made of the same material and utilize the same manufacturing processes.

The sterilized version of the device is a folded gown, packaged with no hand towel, and placed inside a peel pouch. The ENMED™ reinforced surgical gown is also sold as bulk non-sterile, single use

Storage and Disposal:

- Store in dry, clean conditions in original packaging, away from direct sunlight, sources of high temperature, and solvent vapours
- Store within the temperature range 77 to 136° F (25 to 58°C) and with relative humidity below 70%
- Expected shelf life is 2(two) years from date of manufacture when stored as stated (Sterile)
- Expected shelf life is 2(two) years from date of manufacture when stored as stated (Non-Sterile)
- Replace garments if damaged, heavily contaminated or in accordance with local work practice or regulations
- Handle and dispose of contaminated garments with care and in accordance with applicable regulations.

REFERENCE NO	PRODUCT NAME	SIZE	STERILITY
214104001	ENMED REINFORCED SURGICAL GOWN with Towel- Level - 4	S	Non-Sterile
214104002	ENMED REINFORCED SURGICAL GOWN with Towel- Level - 4	M	Non-Sterile
214104003	ENMED REINFORCED SURGICAL GOWN with Towel- Level - 4	L	Non-Sterile
214104004	ENMED REINFORCED SURGICAL GOWN with Towel- Level - 4	XL	Non-Sterile
214104005	ENMED REINFORCED SURGICAL GOWN with Towel- Level - 4	XXL	Non-Sterile
224104001	ENMED REINFORCED SURGICAL GOWN with Towel- Level - 4	S	Sterile
224104002	ENMED REINFORCED SURGICAL GOWN with Towel- Level - 4	M	Sterile
224104003	ENMED REINFORCED SURGICAL GOWN with Towel- Level - 4	L	Sterile
224104004	ENMED REINFORCED SURGICAL GOWN with Towel- Level - 4	XL	Sterile
224104005	ENMED REINFORCED SURGICAL GOWN with Towel- Level - 4	XXL	Sterile

INTENDED USE: The ENMED™ Reinforced surgical gown (ENM-224 sterile and ENM-214 non-sterile) is single use surgical apparel intended to be worn by healthcare professionals to help protect both the patient and the healthcare worker from the transfer of microorganisms, body fluids, and particulate matter.

TECHNICAL D ETAIL

1	COMPOSITION	SSMS REINFORCED TRIPLE-LAYERED MEDICAL PE FILM LAMINATED FABRIC	TOWEL (1 pc)
2	FABRIC WEIGHT	40 GSM	Fabric: %50 viscose, %50 polyester
3	COLOR	BLUE	WHITE
4	SIZE	S-M-L-XL-XXL	40cm x 40cm
5	MADE IN	Turkey	Turkey
6	CARE INSTRUCTION	Single Use	Single Use
7	CLASSES:	Reinforced Gown	Towel
8	SEWING:	Overlock & Ultrasonic	-
9	BELT	Spotwelded	-

Warnings and Limitations


Before use read and understand all user instructions and be sure that the product is suitable for the application and fitted correctly. Product must never be altered or modified. Prior to use inspect the clothing for any damage that could affect its protective function. Replace any damaged clothing. Care should be taken when removing contaminated garments, so as not to contaminate the user with any hazardous substances. If garments are contaminated, then decontamination procedures should be followed prior to removal of the garment. The wearing of chemical protective clothing may cause head stress if appropriate consideration is not given to the workplace environment. Appropriate under garments should be considered to minimize head stress or damage your garment.



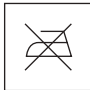
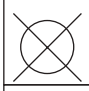





Before using any protective clothing, the wearer must read and understand the user instructions for each product. Specific country legislation must be observed. If in doubt, contact a safety professional. Selection of the most appropriate PPE will depend on the particular situation and should only be made by a competent person knowledgeable of the actual working conditions and the limitations of PPE

Final determination as to the suitability of these products for a particular situation is the employer's responsibility. This information is subject to revision at any time. Always read and follow all User Instructions supplied with your Gown in order to ensure correct operation.

If you have questions contact ENMED Technical Service.

Caution And Precaution



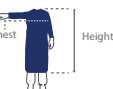
	DO NOT REUSE		DO NOT BLEACH		DO NOT IRON
	DO NOT DRY CLEAN		FLAMMABLE KEEP AWAY FROM SPARKS OF FLAMES		CONSULT INSTRUCTIONS FOR USE
	DO NOT WASH LAUNDERING IMPACTS UPON PROTECTIVE PERFORMANCE		DO NOT MACHINE DRY		CAUTION

Do not use for:

Do not use for chemotherapy drugs, environments with heavy oils, sparks or flame, or combustible liquids; Environments with high mechanical risks (abrasions, tears, cuts); Environments with conditions of excessive heat.

Dimensions:

SIZE	HEIGHT	CHEST
S	47.2 inch	28.9 inch
M	48.4 inch	29.7 inch
L	49.6 inch	30.5 inch
XL	51.1 inch	30.5 inch
2XL	52.3 inch	31.2 inch



LIMITATION OF LIABILITY: Except as provided above, ENMED shall not be liable or responsible for any loss or damage, whether direct, indirect, incidental, special or consequential arising out of the sale, use or misuse of ENMED products, or the user's inability to use such products.

EN ECZA DEPOSU ILAC OZEL SAGLIK HIZ. INS. TAH. A.S.

Address: ATATURK MAH. 31043 SK. KUBAT APT. NO:8 MEZITLI MERSIN , TR 33200

TEL: +90 324 357 56 44 / +90 539 774 03 05

Website/Mail: www.enmedglobal.com / info@enmedglobal.com